


WAX FILTERS / WAX GUARDS



One of the most common problems affecting hearing aids is the ear wax blocking the wax filters. This can result in your hearing feeling dull, or you might even think your aids have stopped working altogether. If ever your aids appear to stop working, try this quick fix first...

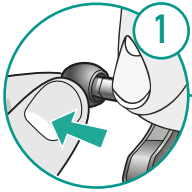



1 Remove the dome by pulling it straight off the receiver. Gently wipe the dome clean with a soft tissue rubbing gently to clean any wax off.

You will see the old white wax filter in the end of the receiver. Take a wax filter stick from the pack. The two ends are different.

A: Filter removal tool (Yellow end) **B:** New wax filter (White end)

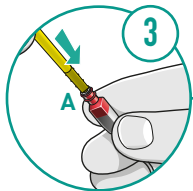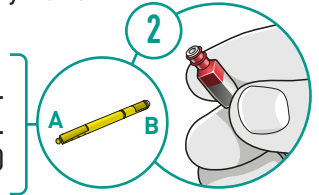


3 Insert removal tool end (A) directly into the existing wax filter and push firmly down keeping the stick straight. Pull the tool straight back out and it will remove the old wax filter

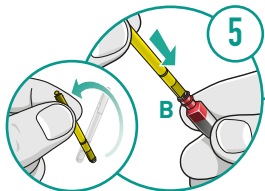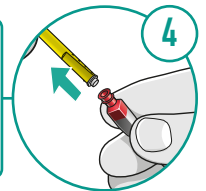

5 Turn the wax stick around and push the new white filter (B) firmly down in to the receiver on your aid. Apply gentle pressure and then remove the stick vertically. The new wax filter should now be in position.

Take the cleaned dome and align the centre of it and push it back onto the receiver. Make sure the dome is securely fitted on by gently pulling it to test it is secure.

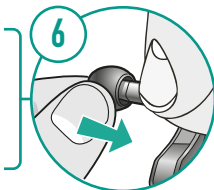

Hopefully, your hearing aids are now back to normal. However, if not, then please call us on **01442 389949** for advice or to book an appointment.

***New wax guards (pack of 8) available at any of our shops for £6.50**

HALO HEARING SOLUTIONS
BERKHAMSTED | BOVINGDON | RICKMANSWORTH

3 Lower Kings Road • Berkhamsted • HP4 2AE | 81b High Street • Bovingdon • HP3 0HP

01442 389949 • info@halohs.co.uk • www.halohs.co.uk