

HEARING AID CHARGING



1

Your new hearing aids will have arrived complete with charging case.

This case should be plugged in at all times and the hearing aids replaced in the case. When fully charged, the hearing aids under normal circumstances will retain enough charge to last you through the day.



2




3

In the evening when you go to bed, simply take your hearing aids out of your ears ensuring that the right hearing aid goes in the right slot and the left into the left. You know they're in the case correctly because the wires will be crossing over each other pointing towards the centre of the case.

If the hearing aids are charging correctly the lights will be flashing green. Once completely charged it will be a solid green however leave them in completely overnight until you're ready in the morning.

4

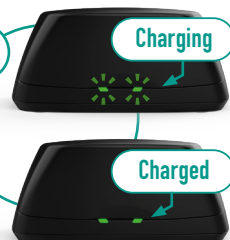

5

When you place them in the case if one of the lights is not flashing or is orange in colour this is signalling a problem with charging.

Remove the hearing aid take a clean tissue and ensure that the contacts at the bottom of the hearing aid are clean and ensure there are no foreign objects in the case itself between the magnetic contact and the charger.

6



Hopefully, your hearing aids are now fully charged. However, if not, then please call us on **01442 389949** for advice or to book an appointment.

 **HEARING SOLUTIONS**
BERKHAMSTED | BOVINGDON | RICKMANSWORTH

3 Lower Kings Road • Berkhamsted • HP4 2AE | 81b High Street • Bovington • HP3 0HP
01442 389949 • info@halohs.co.uk • www.halohs.co.uk